



ANALIZA BIAŁEK



100011235

Laboratory of Analysis

Client: PRIMUSANABOLICS
Order number: 100011235
Name: PRIMUS ANABOLICS PRIMOBOLAN
Number of samples: 1 vial

Batch number: PRIS2352
Recive date: 2025-12-01
Start date: 2025-12-10
End date: 2025-12-12

Sample identification

| Substance | Method of taking a sample | Form |
|-----------------------|-----------------------------|------------------------------------|
| Methenolone Enanthate | Sample provided by customer | Transparent liquid in a glass vial |

Results

| Substance | Sample name | Content | Purity |
|-----------------------|-----------------------------|-------------|--------|
| Methenolone Enanthate | PRIMUS ANABOLICS PRIMOBOLAN | 100,26mg/ml | ----- |

Research comment

The test sample was applied in the amount of 3 μ l to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H₂O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 100,26mg/ml. Height the measurement error is + - 1 mg. The actual state of the sample sent before the test is presented on the photo below

Photo of the test sample



Password

5DL283GF


Analiza Białek sp. z o.o.
KRS: 0000966168
NIP: 8992922970
REGON: 521751875

Headquarters and Laboratory:
ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:
Analysis and research specialist
mgr Arkadiusz Zajęc



Chromatogram

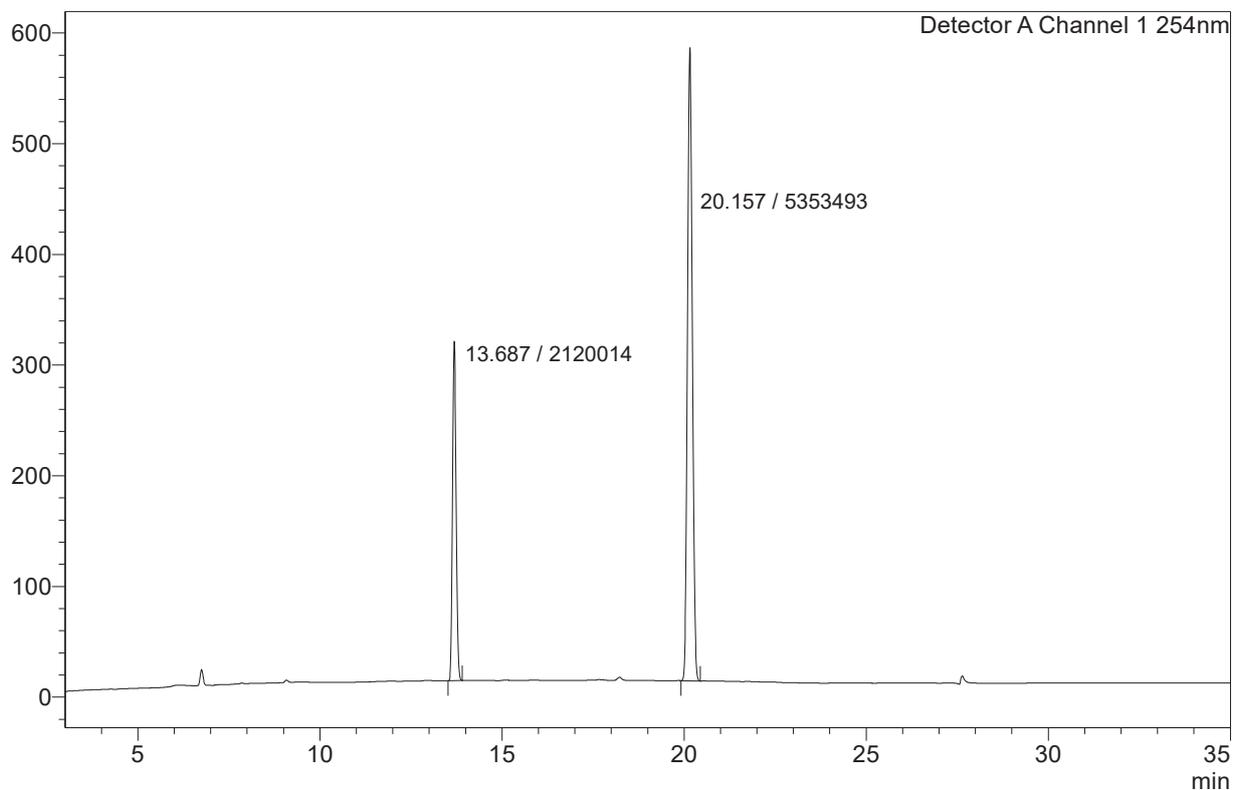

Method

| Time | H2O + 0,1% TFA | ACN + 0,1% TFA |
|-------|----------------|----------------|
| 0min | 60% | 40% |
| 2min | 60% | 40% |
| 27min | 5% | 95% |
| 30min | 5% | 95% |

Password

5DL283GF